DR-HO'S
motionciser
III FIT

Dr. Michael Ho

User manual | P. 2

Brugermanual | P. 30

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Bruksanvisning | P. 37

Manuel d'instructions | P. 16

Ohjekirja | S. 44

Instrukcja | P. 23

Bruksanvisning | S. 51

Thank you

Thank you and congratulations on purchasing your DR-HO'S MotionCiser® Fit, the feel-good massage and passive motion exercise device for your legs and feet while you sit and relax.

We're happy you've chosen to optimize your well-being. Since we want your experiences with the MotionCiser to be fantastic, please read these instructions before use. We wish you years of enjoyment and relaxation!

Including handle to make it easy to move the device about with one hand.

Benefits of DR-HO'S MotionCiser®, from feet to hips:

- helps to loosen foot and ankle joints;
- relax and refresh leg muscles;
- increase circulation in arteries, veins and lymphatics;
- helps to reduce water retention in your legs;
- the unique heat functionality stimulates circulation even more and keep your feet warm and comfortable;
- to help alleviating pain in case of e.g. restless legs.

DR-HO'S MotionCiser® includes:

- 10 levels of intensity for tailor-made effects;
- soothing heat, adjustable from 30°C to 50°C;
- an easy to read LED control panel;
- an easy to use Remote Control, no need to bend over;

- a sturdy and durable design, developed to last;
- silent operation.

DR-HO'S MotionCiser® is recommended for people of all ages

- who lead a sedentary lifestyle due to e.g. study, work or health conditions;
- with (temporary) foot, knee or muscle pain or limited mobility due to surgery, arthritis and such;
- with (chronic) circulation issues, puffiness related to varicose veins e.g.;
- who enjoy a relaxing foot & leg massage in the comfort of their home!

DR-HO'S MotionCiser®:

Low impact exercise with **great** benefits!

PRECAUTIONS AND CONTRAINDICATIONS

Do not use the MotionCiser if you have fracture bones, post-surgery or other health conditions that are not suitable for movement.

Epileptic seizures can be triggered by many factors, including vibration. If you suffer with epilepsy, please consult your doctor before using the MotionCiser.

If you feel any discomfort during or after using the device, stop and rest for a day or two before using it again. If the discomfort continues, stop using the device.

Always consult your doctor with health concerns before use.

IMPORTANT SAFEGUARDS

When using your DR-HO'S MotionCiser®, always observe basic safety precautions.

Warnings

- To avoid risk of electrical shock, never immerse base unit, cord or plug in water/ liquids.
 - Never leave DR-HO'S MotionCiser® unattended while it is in use.
 - This appliance shall not be used or played with by children. Keep the appliance and its cord out of reach of children.
1. Use of the appliance in locations with different electrical specifications may result in damage to the product so please check specifications before operation. Do not use this product with any type of adapter or voltage converter device; this is considered unauthorized modification of the product and as such voids the warranty.
 2. Do not pull, twist, or damage the power cord and don't allow to touch hot surfaces. Periodically inspect cord and plug. Do not operate with a damaged cord or plug. A damaged cord or plug must be replaced by a qualified person in order to avoid a hazard.
 3. If the appliance malfunctions, is dropped or damaged in any manner, contact Customer Service for assistance.
 4. Do not use the device outdoors or in inclement weather, and always operate on a level surface.
 5. This appliance may be used, cleaned and/or maintained by persons with reduced physical, sensory or mental abilities or lack of experience and/or knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved.
 6. Never stand on the appliance but use it for your feet when sitting upright and in a comfortable position.
 7. Switch OFF the MotionCiser and remove the plug from the socket every time you finish using the device and before cleaning it. Proceed as follows: take hold of the plug (not the cable) and pull the plug from the socket. Never pull on the power cable.

Disclaimer: DR-HO's products are not intended to diagnose, treat or cure diseases and disorders. In the event of health problems, consult your doctor for proper diagnosis and treatment.

CONTENTS

- Specifications
- Features
- Use of control panel
- Use of remote control
- Battery change
- Warranty and communications

SPECIFICATIONS & TECHNICAL DATA

1. Control panel
2. Plastic pedal
3. Shell
4. EVA* pedal
5. Socket
6. Remote control
7. Adaptor
8. Screwdriver
9. Connection point for resistance bands (10)

*EVA stands for Ethylene Vinyl Acetate. This is a type of black foam mat used on the two footrest areas. It's very durable and provides a comfortable cushioning effect.

Power.....	60W
Amplitude	4mm
Motor speed	1600r/min – 3500r/min
Level	0-10
Heating.....	30°C – 50°C
Product size	33 x 38 x 12,5 cm
Packing Size	40,5 x 38,5 x 17,5 cm
N.W. / G.W.	3kg / 4,5kg
Fitness standard and class	EN ISO 20957-1 ; HC
Max user weight (sitting position)	40kg


External charger ecodesign information

Manufacturer	Xiamen Eahunt Electronics Co., Ltd. Room 755, Venture Building, Venture Park, Torch High-tech Zone, Xiamen, Fujian, P.R. China
Model: EA048-2400-2000WG Model	EA048-2400-2000WG
Input voltage and AC frequency	100-240V ; 50/60 Hz
Output voltage, current and power	24.0 V ; 2.0 A ; 50,4 W
Average active efficiency	84%
Efficiency at low load (10 %).....	79%
No-load power consumption.....	0,75W

FEATURES

Control panel

1. Increase temperature
2. Reduce temperature
3. Temperature display
4. Power
5. Speed display
6. Speed up
7. Speed down

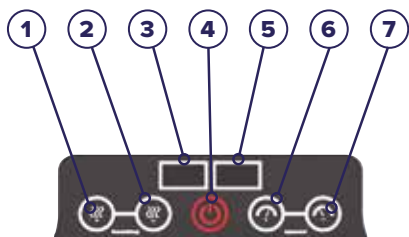Remote control

1. Power
2. Speed up
3. Increase temperature
4. Speed down
5. Reduce temperature

USE OF CONTROL PANEL

- A. Press the  button to turn on, press it again to turn off.

The appliance will cause your foot and ankle to pivot up and down to simulate the motions of walking.

The rapid repeating pivotal motions will help to loosen the joints in your foot and ankle while causing movement in the calf muscles which helps to promote circulation in your legs and feet. The pivotal motion created by the MotionCiser simulates

the motion of the foot, ankle, and calf muscles just like when a person is walking.

The MotionCiser has a preprogrammed mode which automatically turns off after 30 minutes.

Advice:

It is advised to use the MotionCiser for 20 to 30 minutes each time, two to three times a day if you don't go out for a daily walk.

Adjust the speed (level 1 - level 10)

At low speed the MotionCiser will simulate the motions of around 400 steps, at mid speed around 600 steps and at full speed around 1,000 steps. The higher the speed settings, the rapid passive motions will also provide vibrational massage effects.

- B. Press the  button to speed up.

- C. Press the  button to slow down.

Adjust the temperature (30°C, 35°C, 40°C, 45°C, 50°C)

The unique and soothing heat function stimulates circulation even more and also keep your feet warm and comfortable.

- D. Press the  button to set the temperature higher.

- E. Press the  button to set the temperature lower.

USE OF REMOTE CONTROL

- A. Press the  button to turn on, press it again to turn off.

Adjust the speed (level 1 - level 10)

- B. Press the  button to speed up.

- C. Press the  button to slow down.

Adjust the temperature (30°C, 35°C, 40°C, 45°C, 50°C)

- D. Press the  button to set the temperature higher.

- E. Press the  button to set the temperature lower.

BATTERY CHANGE


1. Pull out the battery

2. Change the battery

3. Place the battery back

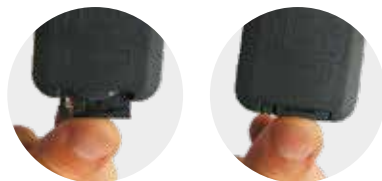WARRANTY AND COMMUNICATIONS

Warranty: one year against factory and material defects. The warranty does not cover faults caused by incorrect handling, damage, assembly or attempted repairs. The same applies to normal wear and tear. This warranty does not affect applicable legal rights.

If you have any comments, questions or concerns, please contact your supplier.

Happy with your DR-HO'S MotionCiser®? Share your posts, tips & tricks! Like & follow us @drhonow & tag your posts

[#MyDrHo](#) [#DrHoLovers](#) [#DrHoNow](#)

SAVE THESE INSTRUCTIONS FOR HOUSEHOLD USE ONLY

THIS MARKING INDICATES THAT THIS PRODUCT SHOULD NOT BE DISPOSED WITH OTHER HOUSEHOLD WASTES THROUGHOUT THE EU. TO PREVENT POSSIBLE HARM TO THE ENVIRONMENT OR HUMAN HEALTH FROM UNCONTROLLED WASTE DISPOSAL, RECYCLE IT RESPONSIBLY TO PROMOTE THE SUSTAINABLE REUSE OF MATERIAL RESOURCES. TO RETURN YOUR USED DEVICE. PLEASE USE THE RETURN AND COLLECTION SYSTEMS OR CONTACT THE RETAILER WHERE THE PRODUCT WAS PURCHASED. THEY WILL TAKE CARE FOR ENVIRONMENTAL SAFE RECYCLING.

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